



Universal Plus®

Electrosurgical and high volume suction/irrigation systems for surgical endoscopy.

All-In-One Instrument System

- Cut
- Cauterize
- Dissect
- Laser
- High Volume Suction
- High Volume Irrigation
- Specimen Retrieval
- Smoke Evacuation



Exclusively Distributed by:



Lap Tech Medical Sdn Bhd

Tel: +603 - 7880 7179

Website: www.laptechmedical.com.my

Email: admin@laptechmedical.com.my

CLEARVIEW™
SMOKE EVACUATION SYSTEM



AIRSEAL™ SYSTEM
INSUFFLATION MANAGEMENT



HELIXAR™ ESU
WITH ARGON BEAM COAGULATION



To learn more about these and other innovative products, call 800-448-6506 or visit CONMED.com

STRAIGHT HANDLE
Electrosurgical Control Options



Hand Controlled Cut and Coag via E/S Buttons with tactile feel and click Hand Controlled S/I



Hand Controlled Cut and Coag Rocker switch Hand Controlled S/I



Foot Controlled Cut and Coag with Universal Cord Hand Controlled S/I

PISTOL GRIP
Electrosurgical Control Options



Hand Controlled Cut and Coag via buttons with tactile feel and click Hand Controlled S/I




Foot Controlled Cut and Coag with Universal Cord Hand Controlled S/I

Universal Plus®

Straight Handle

A small size pencil grip configuration with color coded E/S buttons or rocker switch activation for precise surgeon control. Familiar trumpet valve S/I buttons are concave/convex and color coded for easy identification. The tethered S/I tubing and E/S cord eliminate entanglements and are directed downward at a 45° angle, to minimize surgeon fatigue.

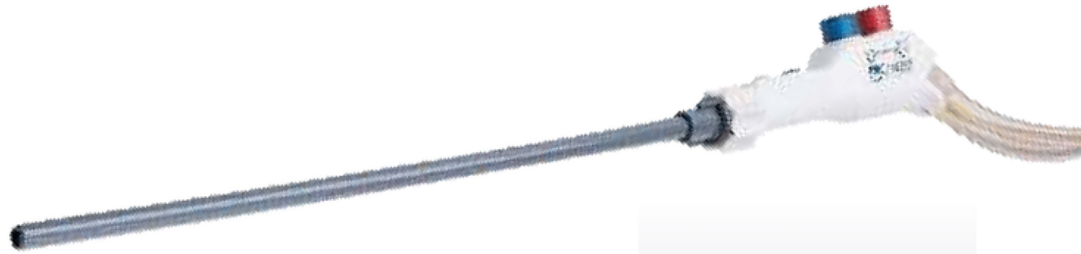

Interchangeable Electrode and Cannula Connection

Pistol Grip


Familiar pistol grip handle design...for surgeon comfort and control. All activation buttons are color coded. The raised S/I controls offer proportional flow control and the convenient cut/coag buttons feature the desired tactile "click" and S/I tubing and E/S power cord are tethered together to prevent entanglement.



Exclusively Distributed by:

Lap Tech Medical Sdn Bhd
Tel: +603-7880 7179

Website: www.laptechmedical.com.my
Email: admin@laptechmedical.com.my

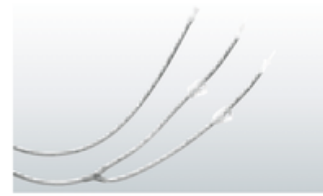

Unique Extendable Lockable Tip Protector Feature

Allows varied tip exposure during electrosurgery. During S/I, it protects viscera while suctioning and directs flow while irrigating.

“Dual-Spike” and Standard Suction Connections

All Universal Plus® Control Handles feature the standard push-on connector for suction canister attachment, and “dual-spike” irrigation line for quick bag changes or dual pump set-ups.

5mm Electrosurgery Electrodes with Suction/Irrigation Optional Stealth ER Eschar Resistant Coating

Available in 27, 32 and 44cm lengths.



Spatula
60-5272-027, 27cm
60-5272-032, 32cm



L-Hook
60-5272-127, 27cm
60-5272-132, 32cm



J-Hook
60-5272-227, 27cm
60-5272-232, 32cm



Needle
60-5272-927, 27cm
60-5272-932, 32cm



Spatula
60-5274-027, 27cm
60-5274-032, 32cm
60-5274-044, 44cm



L-Hook
60-5274-127, 27cm
60-5274-132, 32cm
60-5274-144, 44cm



J-Hook
60-5274-227, 27cm
60-5274-232, 32cm
60-5274-244, 44cm



Needle
60-5274-932, 32cm
60-5274-944, 44cm

Extended insulation at the base of each electrode tip provides greater patient safety, precise surgeon control and reduces possible collateral desiccation. Extendable, lockable tip protector allows coverage of electrode, protecting viscera while suctioning.

Interchangeable Reusable Stainless Steel Suction/Irrigation Cannulas



10/10mm Sump Cannula
60-6002-301, with side holes, 32cm
60-6002-307, with side holes, 15cm



5/5mm Sump Cannula
60-6002-302, 32cm



60-6002-305, no side holes, 15cm

Interchangeable Disposable Suction/ Irrigation Cannulas – Single Use



10/10mm Sump Cannula
60-6001-110, 32cm
60-6001-106, no sideholes, 15cm
60-6001-107, with sideholes, 15cm



5/5mm Sump Cannula
60-6001-203, 32cm
60-6001-205, 44cm



10/10mm “Soft-Tip” Sump Cannula
60-6000-130, 32cm



5/5mm Laser Access Cannula
60-6001-360, 32cm
60-6001-361, 44cm



10/5mm Sump Cannula
60-6001-108, 32cm



Exclusively Distributed by:

Lap Tech Medical Sdn Bhd
Tel: +603-7880 7179

Website: www.laptechmedical.com.my
Email: admin@laptechmedical.com.my

Ordering Information

Description	Catalog Number	Description	Catalog Number
Control Handles		Electrodes with S/I Channel, Extendable, Lockable Tip Protector and Extended Insulation Optional Stealth ER Eschar Resistant Coating	
Pistol Grip, Foot Control E/S, Hand Control S/I, w/10' HiFlo tubing, irrig. bag spike and universal cord, 1/pkg 10/cs	60-6010-000	Spatula, 27cm, 5mm, 1/pkg 5/cs	60-5272-027
Pistol Grip, Hand Control E/S and S/I, w/10' HiFlo tubing, irrig. bag spike, 1/pkg 10/cs	60-6010-002	Spatula, 32cm, 5mm, 1/pkg 5/cs	60-5272-032
Pistol Grip, Foot Control E/S, Hand Control S/I, w/10' HiFlo tubing, "Y" Connector for mechanical pumps and universal cord, 1/pkg 10/cs	60-6010-020	Spatula w/Stealth ER, 27cm, 5mm, 1/pkg 5/cs	60-5274-027
Pistol Grip, Hand Control E/S and S/I, w/10' HiFlo tubing Irrigation "Y" Connector for mechanical pumps, 1/pkg 10/cs	60-6010-022	Spatula w/Stealth ER, 32cm, 5mm, 1/pkg 5/cs	60-5274-032
Straight Handle, Foot Control E/S and S/I, w/10' HiFlo tubing, irrig. bag spike and universal cord, 1/pkg 10/cs	60-6010-001	Spatula w/Stealth ER, 44cm, 5mm, 1/pkg 5/cs	60-5274-044
Straight Handle, Hand Control E/S (Buttons) and S/I, w/10' HiFlo tubing and irrig. bag spike, 1/pkg 10/cs	60-6010-003	"L" Hook, 27cm, 5mm, 1/pkg 5/cs	60-5272-127
Straight Handle, Hand Control E/S (Rocker switch) and S/I, w/10' HiFlo tubing and irrig. bag spike, 1/pkg 10/cs	60-6010-005	"L" Hook, 32cm, 5mm, 1/pkg 5/cs	60-5272-132
Cannulas-Disposable		"L" Hook w/Stealth ER, 27cm, 5mm, 1/pkg 5/cs	60-5274-127
5mm, Sump, 32cm, Single-Use, 1/pkg 5/cs	60-6001-203	"L" Hook w/Stealth ER, 32cm, 5mm, 1/pkg 5/cs	60-5274-132
5mm, Sump, 44cm, Single-Use, 1/pkg 5/cs	60-6001-205	"L" Hook w/Stealth ER, 44cm, 5mm, 1/pkg 5/cs	60-5274-144
10/5mm Sump, 32cm, Single-Use, 1/pkg 5/cs	60-6001-108	"J" Hook, 27cm, 5mm, 1/pkg 5/cs	60-5272-227
10mm "Soft Tip" sump, 32cm, Single-Use, 1/pkg 5/cs	60-6000-130	"J" Hook, 32cm, 5mm, 1/pkg 5/cs	60-5272-232
10mm Sump, 32cm, Single-Use, 1/pkg 5/cs	60-6001-110	"J" Hook w/Stealth ER, 27cm, 5mm, 1/pkg 5/cs	60-5274-227
10mm Sump, no side holes, 15cm, Single-Use, 1/pkg 5/cs	60-6001-106	"J" Hook w/Stealth ER, 32cm, 5mm, 1/pkg 5/cs	60-5274-232
10mm Sump, w/side holes, 32cm, Single-Use, 1/pkg 5/cs	60-6001-107	"J" Hook w/Stealth ER, 44cm, 5mm, 1/pkg 5/cs	60-5274-244
5mm Laser Cannula, 32cm, Single-Use, 1/pkg 5/cs	60-6001-360	Needle, 27cm, 5mm, 1/pkg 5/cs	60-5272-927
5mm Laser Cannula, 44cm, Single-Use, 1/pkg 5/cs	60-6001-361	Needle, 32cm, 5mm, 1/pkg 5/cs	60-5272-932
Cannulas-Reusable Stainless Steel		Needle w/Stealth ER, 32cm, 5mm, 1/pkg 5/cs	60-5274-932
5mm, 32cm, reusable, 1/pkg 1/cs	60-6002-302	Needle w/Stealth ER, 44cm, 5mm, 1/pkg 5/cs	60-5274-944
10mm Sump, no side holes 15cm, reusable, 1/pkg 1/cs	60-6002-305		
10mm Sump, w/side holes 15cm, reusable, 1/pkg 1/cs	60-6002-307		
10mm Sump, w/side holes 32cm, reusable, 1/pkg 1/cs	60-6002-301		
Universal 10' Monopolar E/S cord, Single-Use, 1/pkg 50/cs	60-5168-005		
"Y" Connector set for mechanical irrig. pumps, 1/pkg 10/cs	60-5177-003		
"Y" Set w/dual spike, 1/pkg 10/cs	87101CS		

Exclusively Distributed by:



Lap Tech Medical Sdn Bhd
 Tel: +603 - 7880 7179
 Website: www.laptechmedical.com.my
 Email: admin@laptechmedical.com.my